


Use sectional matrix system



1-1.5mm distance from cusp peaks

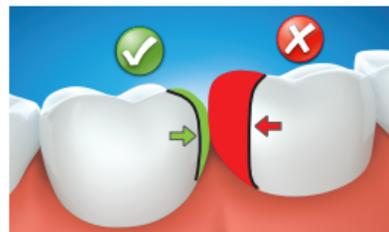

Do not over-contour



1
OPTIONAL STEP:
Apply CAVITY
CONDITIONER (GC) for
10 sec.



2
Rinse and gently dry, do
not desiccate.



3
Shake or tap. Depress
plunger.



4
Mix for 10 sec. Working
time is 1 min. 15 sec.
from start of mix.



5
Insert into Capsule Applicator.
Click twice to prime
capsule.



6
Dispense within
10 sec.



7
Pack and contour.
Avoid moisture
contamination and
desiccation.



8
Final finishing after
2 min. 30 sec.
from start of mix.



9
Finish the restoration
by applying
EQUIA Forte® Coat.
DO NOT AIR BLOW.



10
Light cure for 20 sec.

