

THERMAL PLASTICS

PRO-FORM GLITTER GUARDS

The Pro-form athletic mouthguard line got a little more flash with Glitter Guards. The Glitter Guard material adds a bit of sparkle to the Pro-form lauded line, which delivers protection, affordability and style. Available in packages of 6.



CE

#7965020	Orange	#7965060	Pink
#7965010	Blue	#7965040	Green
#7965030	Yellow	#7965070	Assorted
#7965050	Purple		

PRO-FORM TIE-DYE MOUTHGUARDS

The Pro-form line now boasts Tie Dye designs to the extensive line. The Tie-Dye Laminates are made the same way all Pro-form laminates to ensure the long-proven quality and protection athletes of all levels need. Available in packages of 6.



CE

#7911111	Rainbow	#7965000	Shamrock
#7954000	Patriot	#7934000	Flamingo
#7924000	Wilderness	#7974000	All in
#7944000	Midnight		

PRO-FORM FUN MOUTHGUARDS

The Pro-form line now carries Fun Guard designs to add even more options to the already extensive line of styles and colors. The Fun Guard laminates are made the same way all Pro-form laminates to ensure the long-proven quality and protection athletes of all levels need. The laminates come in packages of 12.



CE

#9600006	Diagonal (red/green/blue/yellow)
#9600929	Polka Dots (red background)
#9598079	Polka Dots (blue background)
#9601039	Polka Dots (white background)
#9599075	The Fangs (black with white fangs)
#9599077	Dripping Blood (black with red drops)

PRO-FORM SAMPLE PACK

Keystone Industries has put together a sample of various thermoplastics. The package contains all of the thermoplastics presently used in the dental industry and comes with technical booklet that details better-known procedures. It contains all of the Pro-form laminated materials, base plates, custom trays, EVA for bleaching and model duplication, retainers, and much more. Each material is labeled with the name of the product and the various applications for which the material can be used. Can easily be cross-referenced to this catalog and the technical booklet.



CE

#7001211	Pro-form Sample Pack
----------	----------------------

PRO-FORM MOUTHGUARD LAMINATES

Pro-form manufactures the only pre-made laminated mouthguard material in the industry. This material has long been considered the leader in custom-made mouthguards. By laminating two layers of soft EVA together, the laminate has excellent tensile strength. All of these features make the Pro-form mouthguard laminate the material of choice for custom-made mouthguards. Square cut to 5"x 5". Also available cut to round in both 120 mm and 125 mm.



1 COLOR SOLID

#9596920	Clear	#9598730	Black
#9600850	Orange	#9600910	Red
#9601040	Yellow	#9598060	Blue
#9600820	Maroon	#9600870	Purple
#9601020	White	#9600180	1 Color, Assorted
#9600800	Green	#9600889	Breast Cancer Pink
#9600886	Pink		

2 COLOR

#9598360	Blue/Green	#9598300	Maroon/Yellow	#9598190	Green/Yellow
#9598240	Black/White	#9598150	Blue/Yellow	#9598320	Black/Orange
#9598210	Blue/Red	#9598130	Red/White	#9598450	Purple/White
#9598110	Black/Yellow	#9598280	Orange/White	#9598380	Green/White
#9598090	Black/Red	#9598260	Blue/White	#9598420	Red/Yellow
#9598340	Purple/Yellow	#9598400	Orange/Blue	#9598710	2 Color, Assorted

3 COLOR

#9598840	Red/White/Blue
#9598900	Canada Red/White/Red
#9600170	3 Color, Assorted

PRO-FORM MOUTHGUARD RESIN SHEETS

The .160 thick Pro-form Mouthguard Resin Sheets are soft for energy absorption, adding more protection. The resin is perfect for mouthguards and niteguards. Available in sheets of 5" x 5" (12.7 x 12.7 cm).



	25/Pkg	300/Pkg		25/Pkg	300/Pkg
Clear	#9597940	#9598010	Yellow	#9597760	#9597930
Red	#9597740	#9597910	Black	#9597660	#9597870
Blue	#9597680	#9597880	Orange	#9597720	#9597900
Green	#9597700	#9597890	White	#9597750	#9597920

PRO-FORM YOUTH MOUTHGUARD

The .120 thick Youth Mouthguard Laminates are great for little kids who are playing sports. Made of soft EVA, for energy absorption and is thin enough to fit the smaller mouth of a child. Available in sheets of 5" x 5" (12.7 x 12.7 cm).



	25/Pkg	300/Pkg		25/Pkg	300/Pkg
Clear	#9611240	#9611280	Green	#9611243	#9611283
Red	#9611241	#9611281	Yellow	#9611244	#9611284
Blue	#9611242	#9611282	Black	#9611245	#9611285

EXTRA-PRO LAMINATES

Keystone Industries now manufactures a triple layered, extra thick (.200) laminated material for mouthguards. Made for optimum protection for athletes like pros, colleges and varsity high school players. Available in solid colors in packages of 12. If you need extra bulk and protection, you might want to try our triple-layered Extra-Pro Laminate.



#9601215	Clear	#9601225	Blue
#9601485	Green	#9601275	Black
#9601515	Red	#9601350	Assorted
#9601535	Yellow		

CE

MOUTHGUARD STRAPS

Keystone Industries offers easy-to-attach helmet straps to be used with Pro-form mouthguards. The strap works great for any sport that has a helmet, such as football, hockey and lacrosse. The strap is especially helpful due to its safety release feature. This allows the mouthguard to break away under less than 12 pounds of pressure.



#9597420	Clear, 25 straps
#9597400	Clear, 12 straps
#9597360	Black, 25 straps
#9597340	Black, 12 straps

BLEACHING LAMINATES

It's become an acceptable dental technique to bleach teeth using carbamide peroxide and a bleaching tray. The Bleaching Laminate is ideally suited for making the bleaching tray and withstanding bleaching products. It is thin, only .020 inches thick once it's formed over the model. It has a foam liner to absorb the bleach, thus decreasing the bleaching time for the patient. The Bleaching Laminate is our best-selling item for bleaching processes.




CE

#9605780 Bleaching Laminates, 50/Pkg

DUAL LAMINATES/E-GASKET

The Pro-form Dual Laminate (E-gasket) is a unique combination of two materials formed into a single laminated sheet. It was made with bruxing cases in mind. The soft liner assures patient comfort and the hard surface assures long-lasting wear. Also, acrylic will bond to the hard surface to allow the dental technician to modify the occlusal surface. Available in 3 different thicknesses: 2, 3 and 5mm.

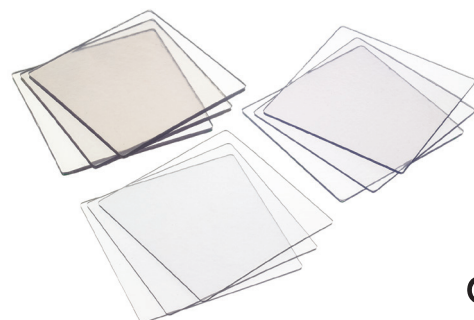

CE

#9604700	2mm E-Gasket, 12 sheets/Pkg
#9603580	3mm Dual Laminates, 12 sheets/Pkg
#9603605	5mm Dual Laminates, 12 sheets/Pkg

NITEGUARDS

This laminate or layered soft/soft material was made specifically for the making of niteguards. The layering gives the thermoplastic a more abrasion-resistant characteristic that will withstand grinding better than traditional soft sheet materials.



CE

#9602550 Niteguards, 12 sheets/box

BASE PLATE MATERIAL

This Base Plate Material is used for denture base or wax bite plates. Made of pure high-impact styrene, the material melts at a low temperature and is easy to work with. Available in 5" x 5" (12.7 x 12.7 cm) sheets and also in round.




CE

	100/Pkg	625/Pkg
0.060 thickness	#9616120	#9616180
0.080 thickness	#9616130	#9616190
0.100 thickness	#9616140	#9616200

COPING MATERIAL

These cloudy, stiff polypropylene sheets are excellent for fabricating temporary crowns, bridges and copings. The material gives good definition and will not adhere to acrylics. Available in 5" x 5" (12.7 x 12.7 cm) sheets.

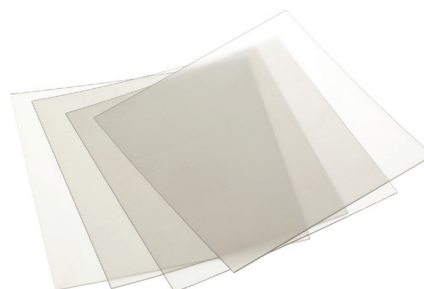

CE

	50/Pkg	625/Pkg
0.020 thickness	#9613470	#9613540
0.030 thickness	#9613480	#9613550
0.040 thickness	#9613510	#9613530

CROWN & BRIDGE MATERIAL

Our Crown & Bridge Material is a thin, clear rigid material that works great for temporary bridges. The material is easy to work with while it does not bond to acrylic. Available in 5" x 5" (12.7 x 12.7 cm) sheets.




CE

#9613680	0.020 thickness, 50/Pkg
#9613700	0.020 thickness, 625/Pkg

RETAINER MATERIAL

This stiff but clear material is ideal for orthodontic retainer. The long-lasting material enables acrylics to stick, but also provides the long-lasting form that retainers need. Available in 5" x 5" (12.7 x 12.7 cm) sheets.

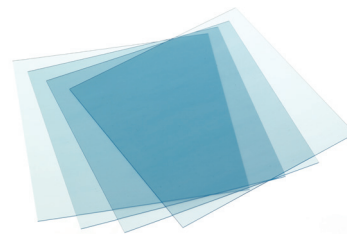

CE

#9614810	Clear 0.030 thickness, 50/Pkg
#9614870	Clear 0.030 thickness, 625/Pkg
#9615025	Super Clear 0.030 thickness, 25/Pkg
#9614816	Super Clear 0.030 thickness, 50/Pkg
#9614875	Super Clear 0.030 thickness, 625/Pkg

PRO-SPLINT MATERIAL

Our new pro split material is a thicker choice for your work in hard and thin splints, surgical trays, retainers, and bleaching trays. It has a less risk of bubbling compared to other splint material due to its thickness, and becomes crystal clear after the heating process is complete.




CE

	25/Pkg	625/Pkg
0.040 thickness	#9640961	#9640965
0.060 thickness	#9640962	#9640966
0.080 thickness	#9640963	#9640967
0.100 thickness	#9640964	#9640968

SPLINT MATERIAL

This stiff clear material is ideal for hard splints, thin splints, surgical trays, orthodontic retainers and bleaching trays. The material sticks to acrylics for various dental appliances. The versatile splint material will surely be used in any lab or office. Available in 5" x 5" (12.7 x 12.7 cm) sheets.

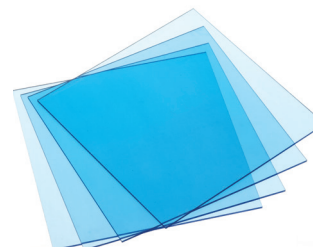

CE

	50/Pkg	625/Pkg
0.020 thickness	#9614800	#9614860
0.040 thickness	#9614830	#9614880
0.060 thickness	#9614780	#9614840
0.080 thickness	#9614790	#9614850
0.100 thickness	#9614820	#9614950

TRAY MATERIAL

The Tray Material is great for custom impression trays. Made of pure high-impact styrene, the material comes in two different thicknesses to provide versatility. Spacer foam is also included with the tray material. Available in 5" x 5" (12.7 x 12.7 cm) sheets and also in round.



CE

#9617320	0.125 thickness, 100/Pkg
#9617350	0.125 thickness, 625/Pkg
#9617330	0.150 thickness, 100/Pkg
#9617360	0.150 thickness, 625/Pkg

SOFT EVA

Keystone Industries' Soft EVA is great for the fabrication of customer fluoride trays, bleaching trays and model duplication. The material is soft, clear and easily forms. Also the material can be trimmed effortlessly. The material is ideal for soft bruxing patients. Available in 5" x 5" (12.7 x 12.7 cm) sheets.




CE

	25/Pkg	300/Pkg
0.040 thickness	#9596980	#9597040
0.060 thickness	#9597090	#9597130
0.080 thickness	#9597150	#9597240
0.120 thickness	#9611240	#9611280

SPACER FOAM

Reusable spacer foam, one size fits all. Available in sheets of 5" x 5" (12.7 x 12.7 cm).

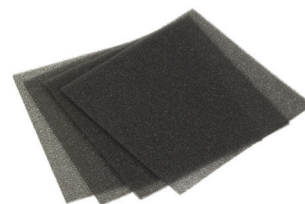

CE

#9617430 Spacer Foam, 25/Pkg

VACUUM FORMERS

MACHINE III

The Machine III has been a staple in Keystone Industries' long-lasting line of highly-recommended vacuum former machines. The price, features and reliability make it the first choice in dental vacuum formers. The machine has more power than any other dental vacuum former, which results in a more accurate and definite dental thermoplastic appliance.

The aluminum construction is long-lasting, eliminating the possibility of rust. The heating element is symmetrically engineered for even heat distribution throughout the thermoplastic material. Made in the USA.

The Machine III is 8-1/4" w x 14" h and it weighs 19 lbs.



#7000330 Machine III, 110 Volt, UL/cUL/CE marked
#7000332 Machine III, 220 Volt, UL/cUL/CE marked

SINGLE CHAMBERED VACUUM FORMER

The Single Chambered Vacuum Former is an economical answer that can still excellently fabricate most thermo-plastic appliances. The all-aluminum construction eliminates rusting and the heating element is symmetrically engineered for uniform heat distribution. The hexagonal post keeps the frame from rotating or sliding. This vacuum former uses 5" x 5" (12.7 x 12.7 cm) resin sheets. Made in the USA.

The Vacuum Former is 9-3/4" x 8-1/4" x 13-1/2" and it weighs 15.8 lbs.



#9594470	110 Volt, UL/cUL marked	#9594380	EU/220 V CE marked
#9594480	220 Volt, CE marked	#9594300	UK/220 V CE marked
	available in non-CE marked	#9594320	Australian/220V CE marked

USA 52 West King Street Myerstown, PA 17067
Phone: (800) 333-3131

Germany Werner-von-Siemens-Str. 14a 78224 Singen
Phone: (+49) 7731 912101 Fax: (+49) 7731 912102

ООО «Денталдепо ЛОТ»
394006, г. Воронеж, ул. Красноармейская, дом 60.
+7 (473) 277-64-77; 271-49-54; 257-25-52
mail@dentald.ru - www.dentald.ru