

The solution for durable bonds



ML Primer

ML Primer is a new, secure, and durable bonding system for veneering composite and all common dental alloys. Based on adhesive monomers, **ML Primer** reliably bonds without requiring extensive processing and is resistant to moisture and thermal strains.

ML Primer is easy to handle and suitable for virtually all indications, as it can be used intra and extra orally.




Proven Products for Better Dentistry

Secure bonding on all dental alloys

ML Primer convinces with an excellent and durable bonding between light or self-curing composites and all common dental alloys.


Especially in the combination with Galvano- and high precious frames

ML Primer provides a secure bonding even after longer time of wearing.

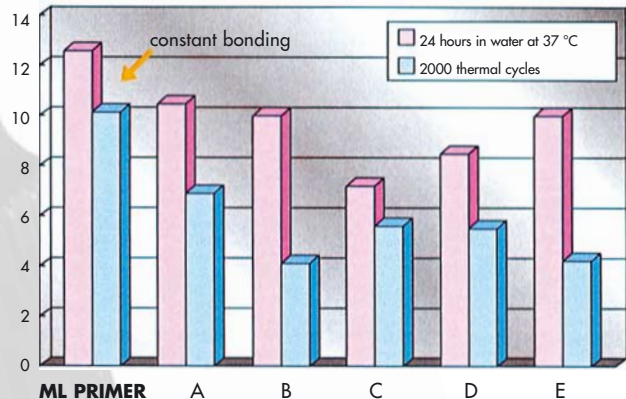


Bonding strength of ML Primer (Mpa)
(high precious alloy in combination with light-curing composite)

Easy and fast application

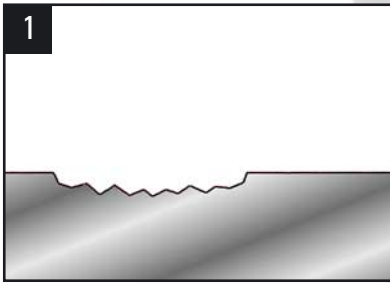

sandblast the surface,
clean and dry

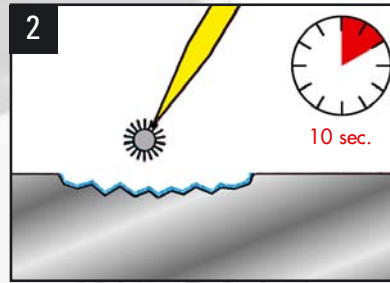

apply ML Primer, leave
for **10 sec.**, ready for
applying composite

The main advantages of ML Primer:

- intra and extra oral use
- unrestricted indications
- easy handling
- fast application
- thin layer
- no additional equipment necessary
- no shadows or dark margins in the bonding area



ML Primer
PN 1740

SHOFU DENTAL GMBH European Head Office

Am Brüll 17 · D-40878 Ratingen · Germany

Phone +49 (0) 21 02 / 86 64 - 0 · Fax +49 (0) 21 02 / 86 64 - 64 · E-mail: info@shofu.de · www.shofu.de

SHOFU DENTAL PRODUCTS LTD. UK Office

Duke's Factory Chiddingstone Causeway · Tonbridge · Kent TN11 8JU · UK · Phone 08 00 / 731 2049

