



initial^{IQ}
LUSTRE PASTES NF

THE FIRST UNIVERSAL 3D PAINTABLE CERAMIC

./'CC'./

INTRODUCING

initial[™] IQ LUSTRE PASTES NF &
initial[™] IQ LUSTRE PASTES NF GUM SHADES





GC INITIAL™ IQ LUSTRE PASTES NF & GC INITIAL™ IQ LUSTRE PASTES NF GUM SHADES:

THE FIRST 3D PAINTABLE CERAMIC TO BRING BEAUTIFUL AESTHETICS TO ALL YOUR RESTORATIVE WORK, JUST BY PAINTING IT ON.

You no longer need complex layering or extra firing cycles to obtain an unsurpassed vitality and a natural glaze for all your restorations. With GC Initial IQ Lustre Pastes NF (based on fine ceramic particles) the steps to complete your crown and bridgework remain identical, regardless of the technique you use.

Do you need to change the color, the brightness, the value, or the surface gloss of your crowns and bridges? You can use the GC Initial IQ Lustre Pastes NF on all ceramics in the GC Initial range, and virtually any ceramic..



Restoration courtesy of Al Hodges CDT, DTG

FEATURES & BENEFITS

GC INITIAL™ IQ LUSTRE PASTES NF

ARE NEWLY DEVELOPED 3-DIMENSIONAL CERAMIC PASTES, WHICH CAN BE APPLIED IN A THICKER LAYER, CREATING COLOR DEPTH AND LIFE-LIKE TRANSLUCENCY.

- With an optimal match to both “low” and “high” CTE range, GC Initial™ IQ Lustre Pastes NF are compatible to almost all types of dental ceramic.
- GC Initial™ IQ Lustre Pastes NF are based on fine ceramic particles allowing a thicker application compared to conventional stains and glazes.
- Due to the unique paste medium, it has a fine thixotropic property for easy and exact applications.
- GC Initial™ IQ Lustre Pastes NF are low fusing, ceramic pastes that do not wear like some stains, glazes or surface sealants.



Restoration courtesy of Bill Marais RDT, DTG

THE DIFFERENCE



WHAT'S SO DIFFERENT ABOUT GC INITIAL™ IQ LUSTRE PASTES?

Most stains and glazes on the market contain oxides that block light



GC Initial™ IQ Lustre Pastes are a unique new ceramic paste that allows light to pass through the restoration

Stains and glazes do not allow you to add morphology, texture, or contacts



GC Initial™ IQ Lustre Pastes allow you to add texture, morphology, and contact by simply controlling the consistency with diluting liquid

Most conventional stains and glazes do not fluoresce at all



Because GC Initial™ IQ Lustre Pastes are indeed a ceramic paste, they fluoresce naturally

Stains and glazes only adhere to the surface of the ceramic and can wear



GC Initial™ IQ Lustre Pastes are low fusing ceramic paste that infuses itself into the restoration

Most stains and glazes must be fired more than once



*GC Initial™ IQ Lustre Pastes are **ONE FIRING!***

Stains and glazes cannot be layered and are not stable



GC Initial™ Lustre Pastes are thixotropic and when left thick can be layered like a ceramic

IMAGINE THE POSSIBILITIES

CAD/CAM

Looking for a way to add extra vitality to your CAD/CAM produced monolithics? GC Initial IQ Lustre Pastes NF is the solution in any case.



Restorations courtesy of Dan Boskocevic

PRESSED CERAMIC



Photoshooting Step by Step: ZTM Michael Brusch – M.B. Dentaltechnik GmbH, Düsseldorf, Germany

Do you prefer to press your ceramics? Start with either the GC Initial IQ - One Body, Press-over-Metal or Press-over-Zircon systems, and then paint the characteristics at a later stage.

This system allows you to characterize your crown and bridgework AND to create a lifelike transition from the crown margin to the gingiva, just by painting.



LAYERED CERAMIC

Achieve the best aesthetic results by using the GC Initial IQ Lustre Pastes NF in the cut-back technique with your pressed or layered work. After painting with GC Initial IQ Lustre Pastes NF, the finishing touches are done with CL-F for high translucency and with the Opal/Incisal masses.



CUT BACK



USE LUSTRE PASTES



OPAL/INCISAL MASSES



FINISHED

*Photoshooting Step by Step: ZTM Michael Brusch –
M.B. Dentaltechnik GmbH, Düsseldorf, Germany*

FULL CONTOUR ZIRCONIA



Restorations courtesy of Masayuki Hoshi

IMAGINE THE POSSIBILITIES

THE FASTEST SOLUTION FOR ALL YOUR MONOLITHIC ZIRCONIUMOXIDE WORK.


- Wide CTE Range in order to work with the complete GC Initial ceramic system and nearly all types of other ceramics or artificial porcelain teeth.
- No complex layering or extra firing cycles required. Fast – Easy – Secure. Ready-to-use pastes with a fine thixotropic property, for easy and precise application.
- Diluting and refreshing liquids to keep the optimal consistency until the end.



*Restorations courtesy of Przemek Seweryniak,
Cosmodent Dental Laboratory, Malmö, Sweden*



CTE

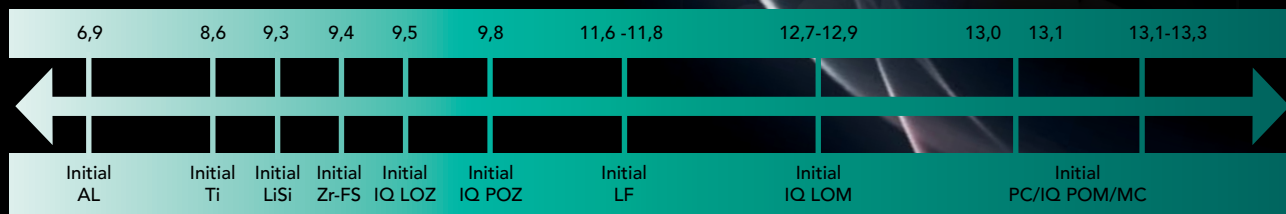

GC INITIAL™ IQ
LUSTRE PASTES NF

GUM SHADES



GC-23
LUSTRE GUM
BASE LIGHT



GC-24
GUM
BASE DARK



GC-34
LUSTRE GUM
INTENSIVE
VIOLET



GC-35
LUSTRE GUM
INTENSIVE
CREAM



GC-36
LUSTRE GUM
INTENSIVE
RED



LP-M1
MODIFIER
WHITE



LP-M2
MODIFIER
RED



LP-M3
MODIFIER
BLUE



LP-M4
MODIFIER
BROWN

ALL YOU NEED TO CREATE A LIFELIKE TRANSITION FROM THE CROWN MARGIN TO THE GINGIVA IN A FEW EASY STEPS.

The five unique gum colored GC Initial IQ Lustre Pastes NF allow an individual adaptation, to match natural gingival shades in natural harmony. You can use the four mixable Gum Modifiers allowing for unlimited creativity and to reproduce the most delicate gum case without any problem.

The GC Initial IQ Lustre Pastes NF Gum Shades are based on a special selection of non-fluorescent, pink colored ceramic particles to perfectly mimic the look of natural gum.



Restorations courtesy of ZTM Michael Brusch – M.B. Dentaltechnik GmbH, Düsseldorf, Germany

GC INITIAL™ IQ LUSTRE PASTES NF V-SHADES

The required V-Shade (value shade) is created by the application of the GC Initial IQ Lustre Pastes NF ranging from "Body A - D". This process, which only takes two minutes, contains everything which is necessary to create V-Shades. Simply use the other available colors to create highly aesthetic results.

With the GC Initial IQ Lustre Pastes NF, you can easily modify your conventional ceramic crown and bridgework. Changing the brightness, value, and color has never been so easy!

The fastest solution for all your monolithic zirconiumoxide work.





BASE



APPLICATION OF BODY SHADES



FINAL RESTORATION OFFERING
COLOR AND GLOSS



BRIGHTNESS



COLOR



VALUE

TIPS & TRICKS

Additional changes in color, transparency and brightness of crown-and-bridge work can be done easily with the GC Initial IQ Lustre Pastes NF. The unique consistency ensures homogeneous modifications and appearance.



BRIGHTNESS

Make your corrections with the effect colors of the GC Initial IQ Lustre Pastes NF Set.

- Increase Brightness > LP Value
- Increase Color Saturations > L A-D
- Increase Gray Value > L 3/4.



BASIS



LP NEUTRAL



COLOR



LP VALUE, L3, LA



GRAY VALUE

Use the GC Initial IQ Lustre Pastes NF to inhibit the depth effects and the transparency of your crown and bridge work build in the layering technique. Sandblast your restoration with 50 micron Aluminum-oxide, damp it and then paint it on with GC Initial IQ Lustre Paste NF.

The requested surface texture/smoothness of the GC Initial IQ Lustre Paste NF layer can be obtained by softly vibrating or condensing the restoration.

Photoshooting Step by Step: ZTM Michael Brusch – M.B. Dentaltechnik GmbH, Düsseldorf, Germany



TESTIMONIALS

"Movies have super heroes to save the day. In real life, we have **GC Initial™ IQ Lustre Pastes** to save our day."

- Step 2thMaster Bay CDT, DTG
CDI SCIENCE, LLC
PINEVILLE, NC



"**GC Initial™ IQ Lustre Pastes NF Gum Shades Set** are a game changer in the way I will handle pink ceramic moving forward, it opens a world of possibilities."

- Bill Marais RDT, DTG
DISA DENTAL STUDIO
PORTLAND, OR



"I was shocked at what I could create with **GC Initial™ IQ Lustre Pastes NF**. This product has made stain and glaze a true art form!"

- Joshua Polansky
NICHE DENTAL STUDIO
CHERRY HILL, NJ




"I use **GC Initial™ IQ Lustre Pastes NF** on nearly every unit that comes out of my lab. It especially makes characterization and glazing of my BruxZir restorations a breeze! The versatility of the paste and natural quality of the colors the kit provides are unrivaled by any other product I have used on the market today. I feel **GC Initial™ IQ Lustre Pastes NF** should be included in every dental artists' arsenal."

- Mike Small
SMALL CERAMICS LLC
LACRESCENT, MN



"I have been using the **GC Initial™ IQ Lustre Pastes NF** kit for a few years now for my CEREC crowns and really like it! The color palette is excellent and gives you a wide variety of options to match teeth. The materials really stay where you place them, so it is easy to characterize with. It has been a very nice addition to my CEREC stain and glaze armamentarium."

- Dr. Rich Rosenblatt
FOREST NORTH DENTAL
LAKE FOREST, IL

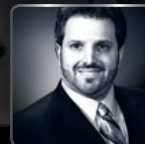

SKUs

GC INITIAL™ IQ LUSTRE PASTES NF

876220 Lustre Paste Diluting Liquid, 8ml
876400 Lustre Paste NF Refresh Liquid, 8ml
876401 Lustre Paste NF Lustre Paste Neutral, 4g
876402 Lustre Paste NF Body Shade A, 4g
876403 Lustre Paste NF Body Shade B, 4g
876404 Lustre Paste NF Body Shade C, 4g
876405 Lustre Paste NF Body Shade D, 4g
876406 Lustre Paste NF Enamel Effect Shade 1 - Vanilla, 4g
876407 Lustre Paste NF Enamel Effect Shade 2 - White, 4g
876408 Lustre Paste NF Enamel Effect Shade 3 - Light Gray, 4g
876409 Lustre Paste NF Enamel Effect Shade 4 - Dark Gray, 4g

876410 Lustre Paste NF Enamel Effect Shade 5 - Light Blue, 4g
876411 Lustre Paste NF Enamel Effect Shade 6 - Dark Blue, 4g
876412 Lustre Paste NF Enamel Effect Shade 7 - INciscio, 4g
876413 Lustre Paste NF Enamel Effect Shade 8 - Olive, 4g
876414 Lustre Paste NF Enamel Effect Shade V - Value, 4g
877043 Lustre Paste NF Glass Spatula
877051 Lustre Paste NF Brush 00
877052 Lustre Paste NF Brush 2
877053 Lustre Paste NF Mixing Dish
877054 Lustre Paste NF Plastic Cover
877078 Lustre Paste NF Set

GC INITIAL™ IQ LUSTRE PASTES NF GUM SHADES

876415 GC Initial IQ One Body Concept LP, NF Gum Shades 4g, G-23
876416 GC Initial IQ One Body Concept LP, NF Gum Shades 4g, G-24
876417 GC Initial IQ One Body Concept LP, NF Gum Shades 4g, G-34
876418 GC Initial IQ One Body Concept LP, NF Gum Shades 4g, G-35
876419 GC Initial IQ One Body Concept LP, NF Gum Shades 4g, G-36
876420 Stain Modifier, for LP NF, 3g LP-M1-White
876421 Stain Modifier, for LP NF, 3g LP-M2-Red
876422 Stain Modifier, for LP NF, 3g LP-M3-Blue

876423 Stain Modifier, for LP NF, 3g LP-M4-Brown
877051 IQ LP NF Brush 00
877052 IQ LP NF Brush 2
877053 IQ LP NF Mixing Dish
877054 IQ LP NF Plastic Cover
876220 IQ LP NF Diluting Liquid, 8mL
876400 IQ LP NF Refresh Liquid, 8mL
877088 GC IQ One Body Concept, LP NF Gum Shade Kit



Restoration courtesy of Bill Marais RDT, DTG



TOUCH*of*GENIUS™

