

GC KALORE MULTI SHADE BUILD-UP GUIDE

CLASS
I & II



CLASS
III



CLASS
IV



JUNIOR



WT

ADULT



DT

NT

SENIOR



GT

CT

CLASS
V



- Translucent (Incisal edge)
- Translucent (Enamel)
- Universal
- Opaque (Opaque Dentin)

Available shades


Opaque



Universal



Translucent

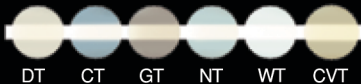

CV: Cervical (B5)

CVD: Cervical Dark (B7)

CVT: Cervical Translucent

DT: Dark Translucent

CT: Clear Translucent

GT: Gray Translucent

NT: Natural Translucent

WT: White Translucent

BW: Bleaching White

XBW: Extra Bleaching White

OBW: Opaque Bleaching White

OXBW: Opaque Extra Bleaching White

Shades by number

	A1	A2	A3	A3.5	A4	B1	B2	B3	C2	C3	D2	CV	CVD	BW	XBW
Opaque	OBW	AO2	AO3	AO3	AO4	OBW	AO2	AO3	AO3	AO4	AO2	AO4	AO4	OBW	OXBW
Universal	A1	A2	A3	A3.5	A4	B1	B2	B3	C2	C3	D2	CV	CVD	BW	XBW
Translucent	WT	WT	DT	DT	DT	WT	WT	DT	DT	DT	WT	DT	DT	WT	WT

GC KALORE

13001248
260392DK